


ELECTRIC WHEELCHAIR



PXJFK005BK

OPERATING MANUAL

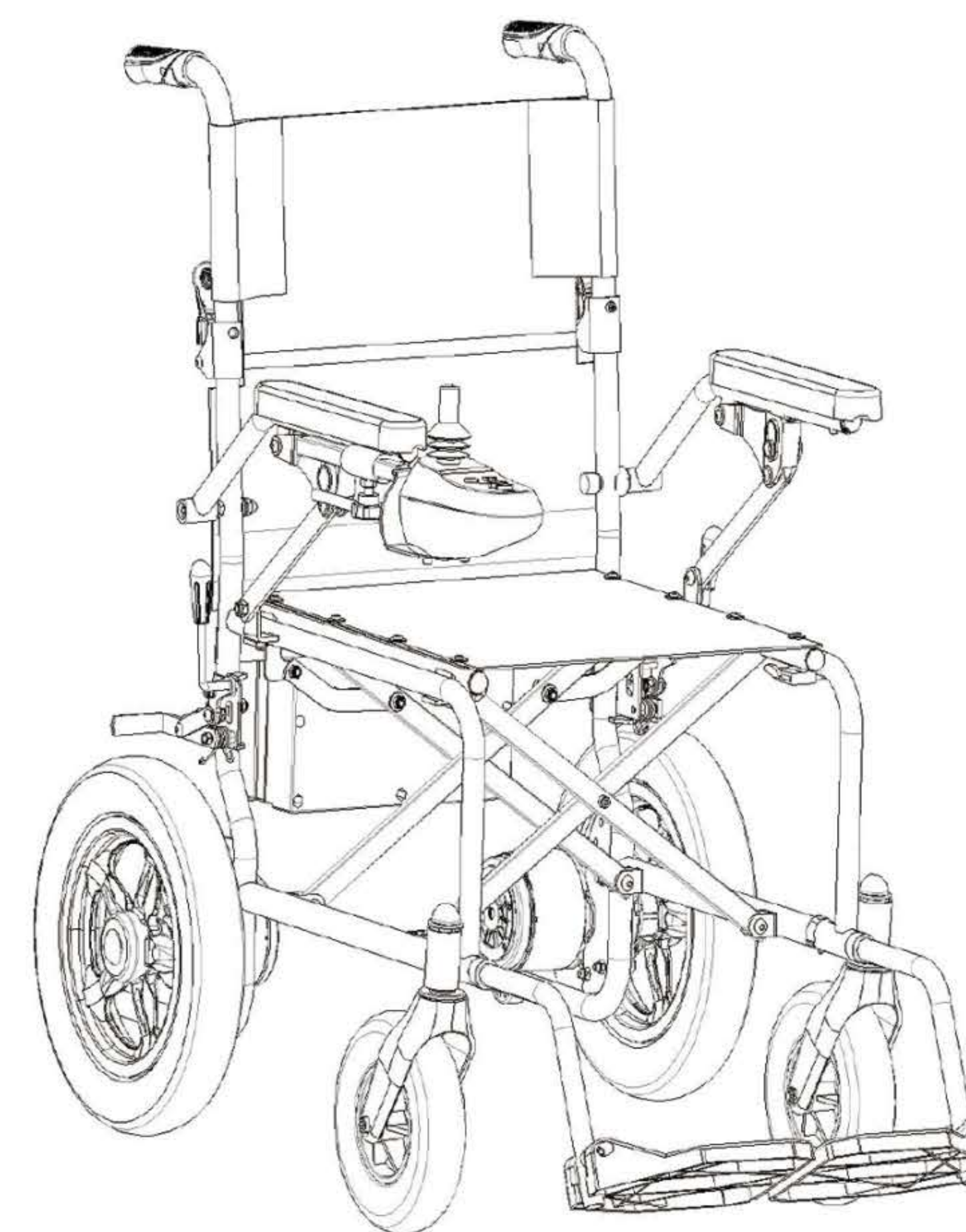

For the safety of you and others, please read the manual carefully!
Buy and use it under the guidance of medical staff!
No driving on the highway!

Catalogue

1. Customer notification
2. Safety Instructions
 - 2.1 Precautions before use
 - 2.2 Notes for escorts
 - 2.3 Safe operation instruction
 - 2.4 Special Instructions on product labels, chargers and use
3. Product overview
4. Installation and use
 - 4.1 Unpacking, assembly and debugging
 - 4.2 Use
 - 4.3 System diagnosis
 - 4.4 Control system
5. Maintenance
 - 5.1 Cleaning and routine inspection
 - 5.2 Regular maintenance records
6. After-sales service
 - 6.1 Product term
 - 6.2 Product warranty period
7. Electromagnetic compatibility
8. Packing list
9. Power supply conditions and environmental requirements
10. Labels, Packing Logo Design

1. Customer notification

Scope of application: The Electric wheelchair is a motor driven, indoor and outdoor transportation vehicle with the intended use to provide mobility to a disabled or elderly person limited to a seated position.

First of all, thank you for choosing to buy and use our electric wheelchair, thank you for your trust and support to our company! Please read our instructions carefully before use!

2. Safety instructions

★ Be sure to read:



Note

Please follow the instructions and specifications in this manual when operating the wheelchair to ensure the safety of the product. Be sure to abide by national and local road traffic regulations and drive carefully.

We only provide the repair and service described in this manual, please use factory parts. Please note that our instruction manual contains many precautions, any vehicle can cause injury if used incorrectly, careless driving may endanger your own safety, but also affect the safety of others. Please follow our guidelines and use your wheelchair properly. When driving outdoors, please obey the traffic rules and follow the points of caution in this manual.

2.1 Precautions before use

2.1.1 Every time you get into or get out of the wheelchair, don't stand on the foot pedal, step on the ground.



Warning

In and out of the wheelchair, do not stand on the pedal!

2.1.2 Practice in a safe and wide area such as a park before you are used to driving.

2.1.3 Make sure the user wears a seat belt to ensure safety.

2.1.4 Make sure to go forward up and down hills. Please proceed at the lowest speed when going down hills.

2.1.5 Please avoid crossing high steps. This product can cross 1.6 inch high steps with

19.7 inch run-up. When crossing the steps, please be sure to walk at a right Angle to the steps.

2.1.6 If you choose to go downhill with your back to the side, please make sure that there is an attendant to take care of you.

2.1.7 Do not use this wheelchair to carry goods and traction purposes. It is prohibited to stand on the seat cushion and jumping.

2.1.8 The circuit design of this wheelchair is based on the safety of the user, and it is strictly prohibited to modify the circuit without permission

 Warning	The user should not drive or use the wheelchair in any motor vehicle because the user may be thrown from the wheelchair and injured in an emergency or emergency braking. When in any motor vehicle, the user shall move from the wheelchair to the normal seat of the motor vehicle and secure the wheelchair in the vehicle.
---	--

2.2 Notes for escorts

2.2.1 Please make sure that the user's feet are securely on the pedals and the clothes will not touch the wheels.

2.2.2 Please keep the controller is turned on with electronic brake when going up or down the hill. Do not turn off the controller.

2.2.3 Please do not brake too hard or too fast in case the user is ejected from the wheelchair.

2.2.4 Make sure the hand brake is on position when the wheelchair is not being driven.

2.2.5 Maximum working load:260lb

 Advice	* Environmental protection: Waste batteries and vehicles should be recycled in accordance with relevant national guidelines. * Speed limit: According to national regulations, the speed of electric wheelchairs should be controlled below 6km/h. Please comply with relevant regulations to protect your own safety.
--	--

2.3 Safe operation instruction

Operation mode of different road surface

Step Up and Down (see Figure 1)

- Slowly approach steps or obstacles and stop at the right Angle;
- Place the front wheel close to an obstacle or step and lean against the step;
- Keep the front wheel at right angle to the obstacle and keep the direction forward;

- Keep the direction and drive slowly forward until the rear wheels are cleared the obstacle and the maximum obstacle height is 4cm.

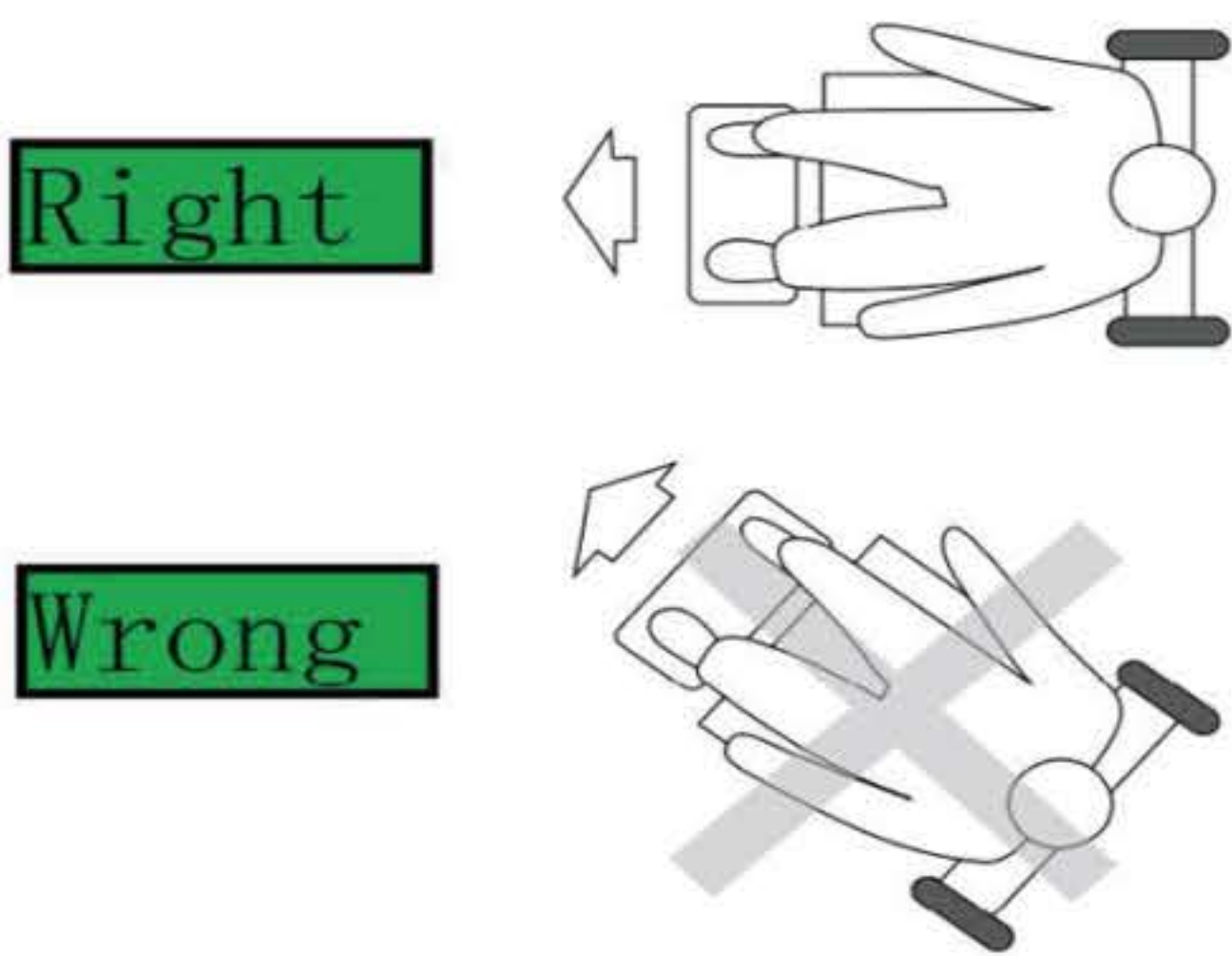


Figure 1

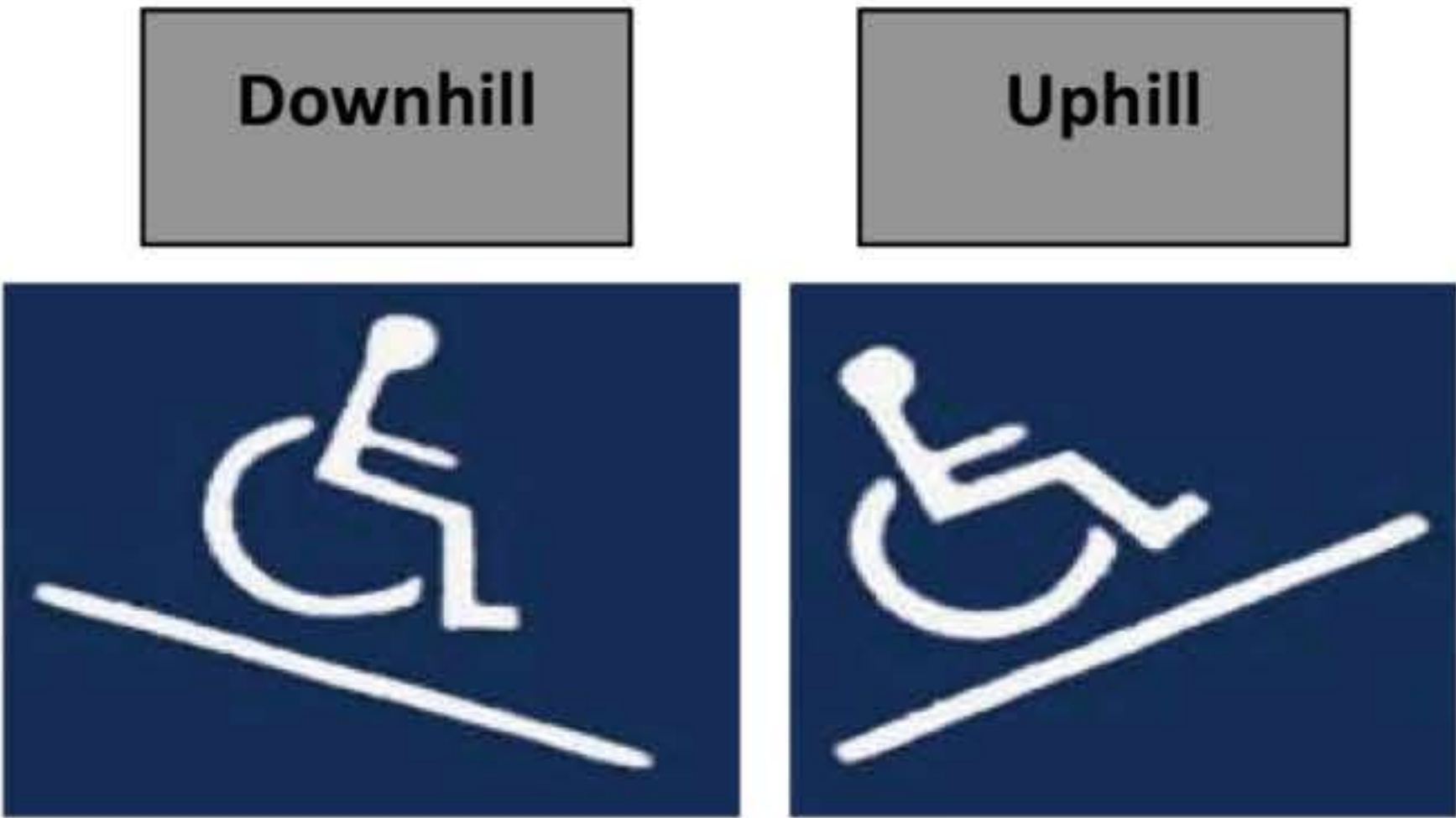


Figure 2

Downhill (see Figure 2)

- The downhill speed should not set to the lowest gear, drive with the lowest speed and keep within a controllable range;
- Do not swerve or brake suddenly;
- Do not try to go downhill on slippery roads;

Uphill

- Keep the back straight when going uphill;
- Do not try to push your torso out of the wheelchair seat range while climbing;
- Drive straight when going uphill. Don't zig-zag or make sudden turns;
- Do not attempt slopes that are too steep and no more than 10 degrees is recommended.

2.4 Special instructions on product labels, chargers and use:

2.4.1 Interpretation of graphics, symbols, abbreviations and other contents used in medical device labels:



Note



Type B application section

2.4.2 The Wheelchair and battery charger comply with GB/T 18029.21;

2.4.3 In addition to the contents specified in ISO 7176-15, the manual contains warnings on the related risks caused by electromagnetic radiation source interference in the use environment of wheelchairs, including but not limited to:
 -- mobile phones;

-- Electronic goods anti-theft system.

3. Product Overview

【Name】 Electric Wheelchair

【Model】 PXJFK005BK

【Scope of application】 The Electric wheelchair is a motor driven, indoor and outdoor transportation vehicle with the intended use to provide mobility to a disabled or elderly person limited to a seated position.

【Structure and composition】 It is composed of frame, leg rest, front wheels, rear wheel, seat and backrest cushion, armrest, controller, motor, battery, charger.

【Product Performance】

Description		Performance index (outdoor)
Max speed, km/h		≤6.0
Running braking performance	Level road braking, m	≤1.5
	Maximum slope safe braking, m	≤3.6 (6°)
Slope parking performance		≥9°
Static stability	Forward leaning	≥9°
	Backward leaning	≥9°
	Side leaning	≥9°
Dynamic stability	Backward dynamic stability on the slope	≥6°
	Backward dynamic stability on the slope	≥6°
	Side leaning dynamic stability on the slope	≥6°
Height of obstacle, inch		≥1.6
Ditch width, inch		3.93


Climbing ability (Angle)	≥6°
Minimum turning radius, inch	≤47.2

【Product Features】

1. Frame:carbon steel; 8inch front wheel; 12inch rear wheel;
2. Motor power: DC24V/250W; dual-motor drive;
3. Rated voltage: D.C.24V;
4. Size(L*W*H): 40.9inch*26.7inch*35inch;
5. Net weight: 65 lb (±10lb);
6. Power supply: 110V/220V 50Hz; Maximum output current:2A
7. Battery: Lead-acid battery 24V12AH;
8. Maximum working load:260 lb

【Name and illustration of wheel chair parts】

- 1.Push handle; 2.Folding lever for the back; 3.Handrail; 4.Up and down button for the handrail; 5.Joystick; 6.Battery; 7.Hand brake; 8.Rear wheel; 9.Motor; 10.Guard wheel; 11.Front wheel; 12.Backrest cushion; 13.Seat cushion; 14.Frame; 15.Middle support frame; 16.Bending tube for Footrest; 17. Footrest belt; 18.Footrest; 19.Charger

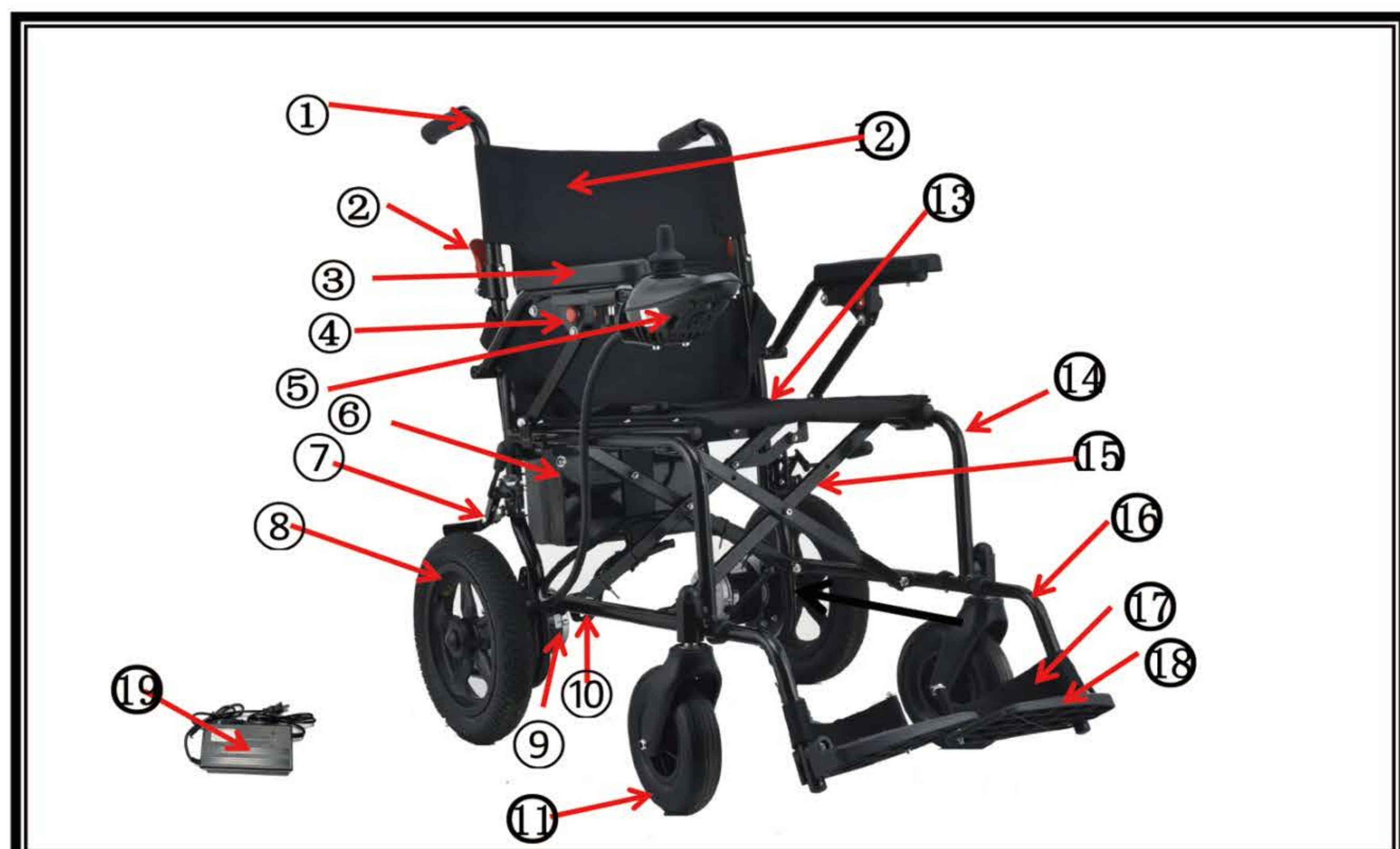

Name, status, and function of each part of the controller:

Controller:

- ① Power: 20% of electric power per grid
- ② Switch: Power on or off
- ③ Horn: Use the horn
- ④ Speed: The current speed
- ⑤ Downshift: Reduce driving speed
- ⑥ Up shift: Increase driving speed
- ⑦ Joystick: Direction, speed and parking



4.Installation and use:

4.1 Unpacking, assembly and debugging

4.1.1 Unpack and check if the wheelchair and accessories are in good condition:

4.1.1.1 Use a small knife to cut open the packing tape and sealing tape. (as shown in Figure 3);



Figure 3

4.1.1.2 Check whether the items in the box are damaged, if there is damage, please take a picture and contact the merchant as soon as possible (as shown in Figure 4);

4.1.1.3 Take out the parts and check if there is any missing or damaged. If there is any missing or damaged, please take photos and contact the merchant as soon as possible [as shown in Figure 5];



Figure 4



Figure 5

4.1.1.4 Two people are required to grasp the steel parts at both ends with both hands, and at the same time, lift them out with force and put it down lightly;

4.1.1.5 Check again whether the wheelchair and its accessories are damaged or missing. If there is any missing or damaged, please take photos and contact the

merchant as soon as possible;

4.1.2 Quick installation and use:

4.1.2.1 Controller installation (for the user or the attendant): Insert one end of the telescopic tube of the controller into the controller fixing tube under the tube of the right handrail, adjust the controller to the appropriate position, and then lock the knob. The power connection wire is fixed backward to reach a straight line [as shown in Figure 6 and Figure 7];



Figure 6



Figure 7

4.1.2.2 Unfold the wheelchair (for the attendant): Press down firmly with both hands in a balanced way until the vehicle is fully open [FIG. 8 and 9]. If still not fully open, press the battery hanger [FIG. 10].



Figure 8



Figure 9



Figure 10

4.1.2.3 Pedal installation (for patients or caregivers): Insert the hanging foot into the frame pipe and thread the hex head screw through the vehicle Frame and footrest, finally position with the hexagonal hole on the frame, and tighten the knob for fixation. As shown in Figures 11 and 12;



Figure 11



Figure 12

4.1.2.4 Battery installation (for the user or the attendant): Hang the battery box on the battery hanging rod of the vehicle, and connect the red plugs tight enough. When charging the battery independently, use both hands to press the red plugs and pull them apart [as shown in FIG. 13, 14, 15, 16].

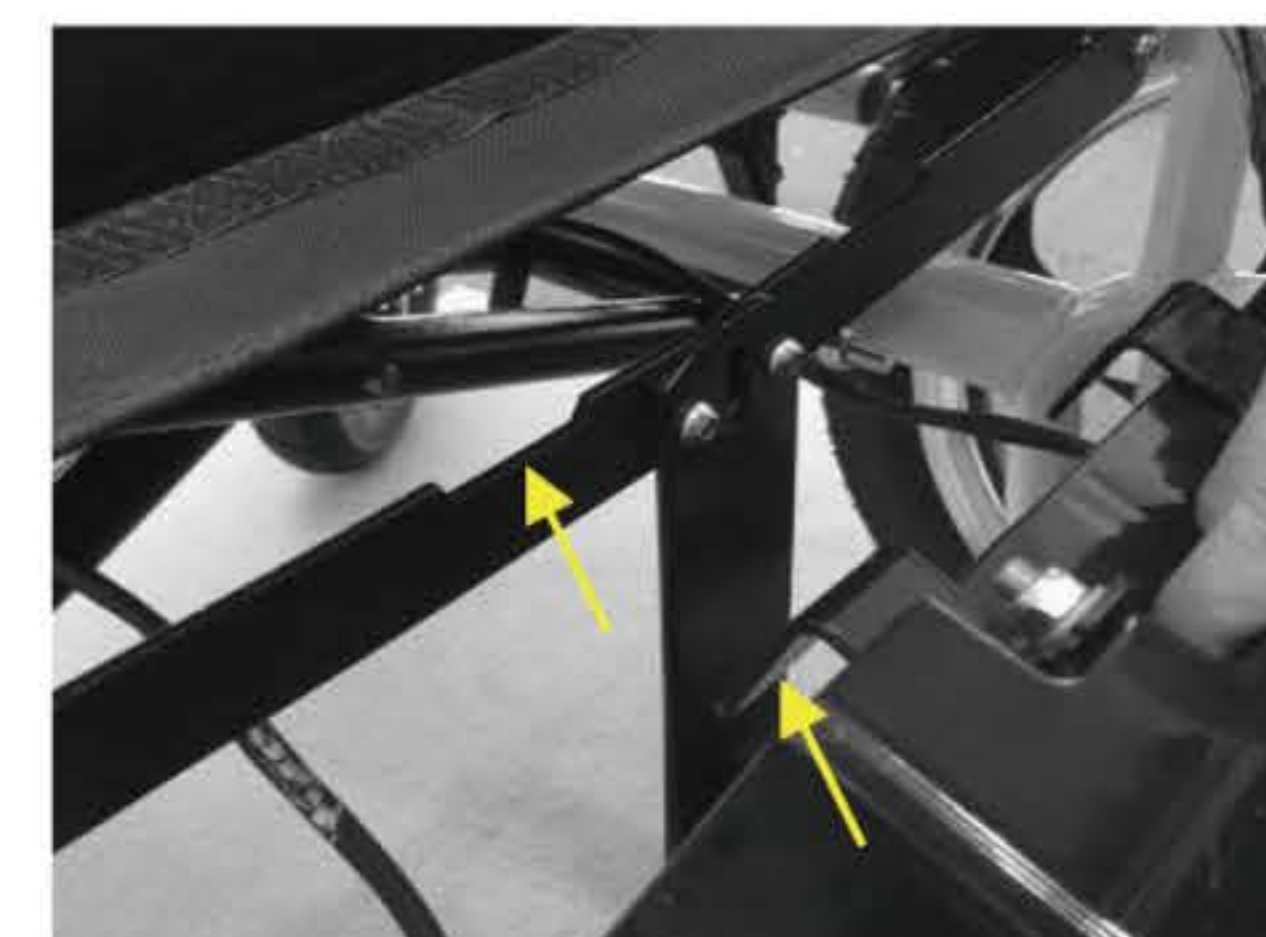


Figure 13

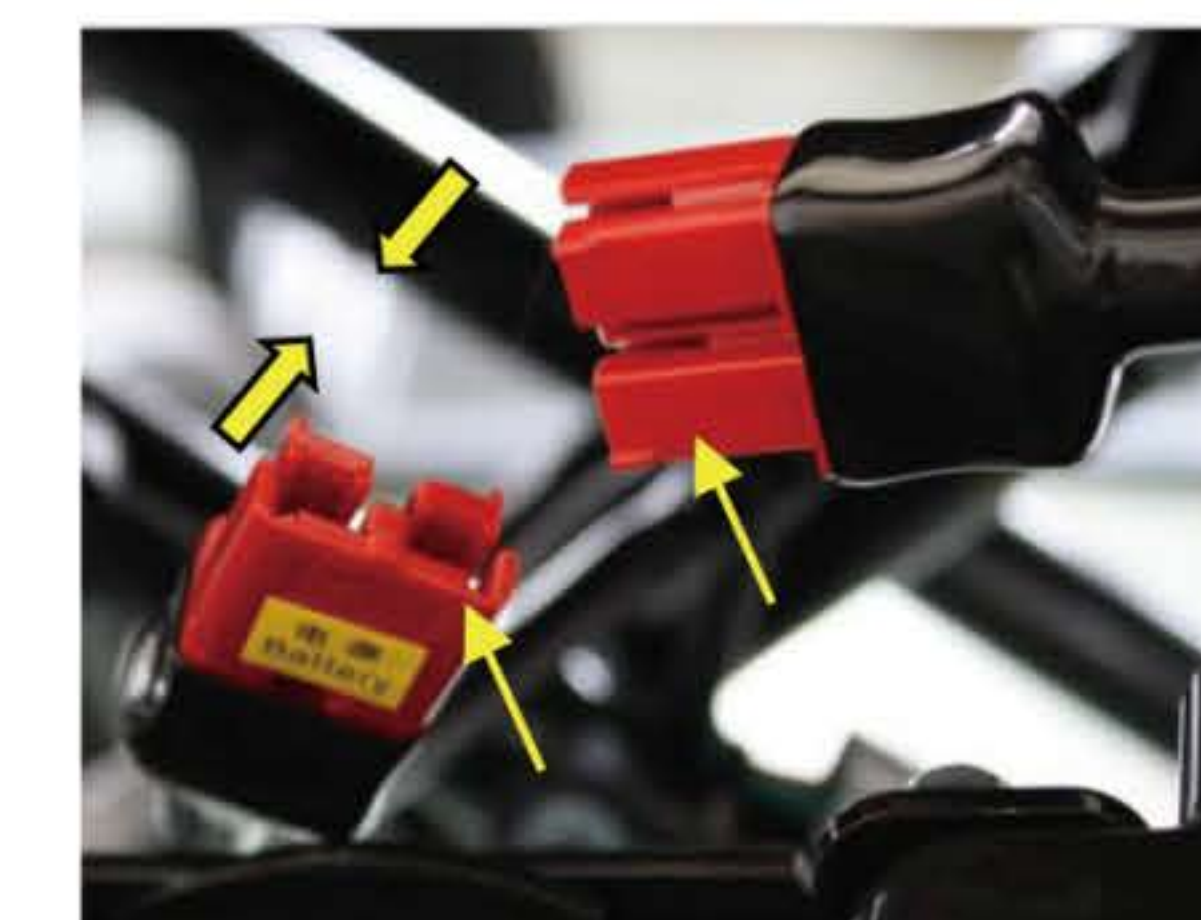


Figure 14



Figure 15



Figure 16

4.2 Use

4.2.1 Backrest adjustment (for the attendant): Hold the push handles of the wheelchair with both hands, lift them upward until the backrest is fully expanded, and make sure that the pin is stuck in the hole of the handle bar with click sound. If you need to fold then press the lever on both sides, please check it in detail to ensure the safety of use [FIG. 17, 18, 19].



Figure 17



Figure 18

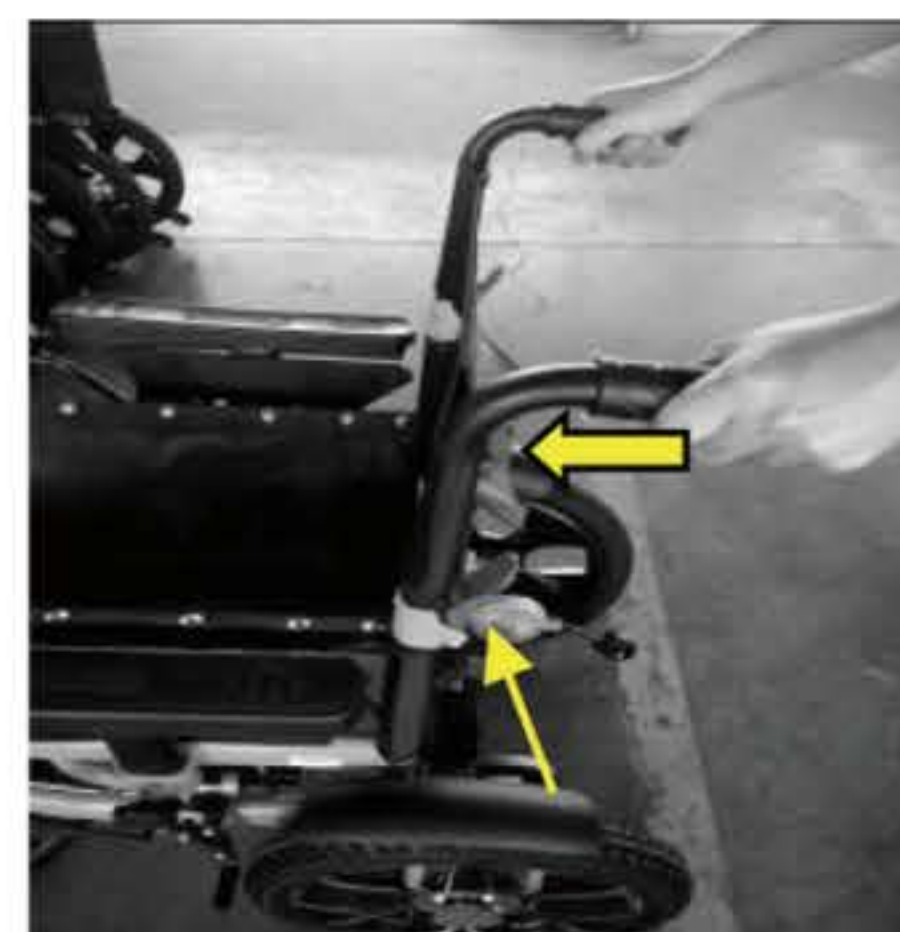


Figure 19

4.2.2 Use of electronic brakes (for the user or the attendant): When the wheelchair needs to be used in "manual mode", it only needs to turn off the controller. In manual mode, the attendant is required to take care of the wheelchair. In case of accidents, the hand brake is in use when the wheel chair is stop. "Electric mode" is started when the controller is power on (if the wheelchair is faulty or the battery power is insufficient or the chair is in manual mode, please manually break the chair)

Manual mode:
The controller is turned off.

Electric mode: The controller is power on, Loose the hand brake, electronic brake is in use.



Note

Please pay attention to the tire condition, when the clutch is disengaged (automatic electric brake disengagement), When switching from electric mode to manual mode, an escort must be present to assist.

4.2.3 Use of the handbrake (for the user or the attendant): hold the handbrake grip and pull back, the brake pad touches the rear tire, fix the handbrake, and the wheel chair stops moving. Please prevent sliding or touching by mistake. The rear wheel is inflatable tire, it is recommended to inflate according to the air pressure marked on the side of the tire, please be careful not to touch the wheel with your fingers, please prevent anyone from putting their hands close to the wheel and causing injury [as shown in Figure 20].



Figure 20

4.2.4 Press the red lock button and pull up the armrest to complete the rear lift of the armrest (as shown in Figures 21 and 22); Press the fixed hand forward until the armrest moves forward until the red lock button clicks and automatically locks.



Figure 21



Figure 22

The round charging connector (Canon head) is directly connected to the controller: charge in the slot below the joystick [as shown in Figure 23] [charger optional]; Please refer to the instructions attached to the charger for further information.



Figure 23

4.2.6 Timing of charging



Advice

For the round charging head charger for the controller, make sure that the controller is turned off for charging, and the actual charging time should not exceed 12 hours;

4.2.6.1 Continuous use for more than one hour, or more than one quarter of the electric quantity is used when the indicator is orange;

4.2.6.2 It is recommended to charge it every day as long as it is used;

4.2.6.3 If the electric wheelchair is not used for a long time, please charge it every 1-2 months;

4.2.6.4 When the remaining power cannot reach the expected distance;

4.2.6.5 Follow the instructions below to complete the charging procedure.

Before charging, please read the instructions provided with the charger and select the

correct voltage:

Step 1: Check that the charger slot is not blocked;

Step 2: Make sure the wheelchair switch is turned off;

Step 3: Connect the output plug of the charger Insert it into the charging slot and make sure it is firmly inserted;

Step 4: Insert the power plug on the charger into the charging slot on the battery or on the controller (see Figure 23), the charger indicator will turn orange (charging), the whole charging process takes about 8~12 hours;

Step 5: When the charging indicator turns green, it indicates that the battery is fully charged;

Step 6: Turn off the power supply connected to the charger, disconnect the output plug from the wheel chair controller charging slot.



Warning

Using it without confirming full charge will reduce battery life. If the accumulated mileage (the mileage between the two charging times) often exceeds half of the maximum endurance, the battery life will be significantly reduced; when charging, personnel should not be on the electric wheelchair. The company will not be responsible for any product failure or damage caused by failure to follow the above battery maintenance conditions or by self-replacement.



Advice

Do not expose the charger to the outdoors or heat sources, such as radiators, fire sources, sunlight.

Do not move the wheelchair until the battery charger plug and the main power plug is disconnected from the wheelchair and the power source. When charging is complete, the charging indicator light will turn into "green light", try not to stop charging before charging is complete.

Unplug the charger after charging, but the longest charging time should not exceed 12 hours, there is a danger of overcharging.

The charging time is related to the external temperature, and a longer charging time is required in winter.

Please observe the following rules to avoid danger when charging:

1. It is strictly forbidden to change or modify the charging circuit without permission. The company will not be responsible for any

	accidents caused to products or personnel; 1. Do not disassemble or modify the charger; 2. The charging place should be well ventilated, do not expose to the sun or humid environment for charging, and do not cover any tarpaulins or objects when charging; 3. When the charger is working, there will be fan noise, please feel free to use it. The function is to dissipate heat, but the temperature of the charger shell will still rise slightly, which is a normal condition; 4. This charger has no waterproof function; 5. Do not place the charger over flammable objects, such as oil, pedals or seats... And so on; 6. Lead-acid batteries have no memory effect. Keeping good charging habits will help prolong battery life. Excessive discharge (less than 1/5 of the residual power) will shorten battery life.
--	---

 Warning	1. Please keep away from the fire source when charging, which may cause the battery to catch fire or explode; 2. Please do not smoke when charging because hydrogen will be generated, please charge in a well-ventilated place; 3. Do not install or remove the charging socket when your hands are wet or the socket is damp. This may cause electric shock.
--	--

4.2.7 Charger

4.2.7.1 Charger parameters: 1.Input voltage: AC110~240V; 2.Output voltage: 24V/DC; Output current: 2A;

4.2.7.2 During charging, the charging indicator will turn orange and when the charging is completed it will turn green;

4.2.7.3 Charger operation is based on the attached instruction;

4.2.8 Battery

4.2.8.1 Do not store the wheelchair or the battery in an environment lower than minus 10 degrees Celsius or higher than 50 degrees Celsius. The above environment will lead to overheating and damage the battery or reduce the battery life;

4.2.8.2 This vehicle uses a maintenance-free battery without replacing or supplementing battery fluid;

4.2.8.3 For other battery-related matters, see the precautions on the battery box sticker.

 Warning	1. Do not open the cover of the battery or get close to the fire, which will cause explosion and danger. 2. Damaged batteries have environmental protection and safety issues, so please recycle them in an environmentally friendly manner or notify the dealer or the company to deal with them.
---	---

4.2.9 Battery cleaning

If the battery is contaminated with dirty water, battery acid or other dust, the battery will discharge rapidly, so please follow the steps below to clean the battery.

4.2.9.1 Turn off the power;

4.2.9.2 Wipe the wheelchair body around the battery with a clean cloth;

4.2.9.3 Wipe the battery with a clean cloth. If there is white powder in the connector, remove it with a copper brush;

 Note	1. Make sure the connector is properly installed; 2. Do not use the battery of this vehicle to supply power to communication equipment or other equipment; 3. The capacity of the battery will change with the external temperature, and the endurance will be shorter in winter.
---	---

4.2.9.4 Battery replacement

 Note	1. For battery replacement, only use the same battery or the same type of battery recommended by the manufacturer. If you need to replace the battery due to special factors, please be sure to read the instructions in this manual. If you have any questions about battery replacement, please contact our company or the dealer you purchased it from; 2. Before replacing or installing the battery, be sure to turn off the power of the controller; 3. Please do not touch any metal objects with the battery connector; 4. Please read the warning label on the battery carefully; 5. Do not dissect the battery to avoid the leakage of the battery fluid and injury to the skin or eyes. If the skin and clothes are stained with the electrolyte, rinse with plenty of water. If it splashes into the eyes, rinse eyes with clean water for at least 15 minutes and seek medical attention in time; 6. For protection, it is recommended to wear rubber gloves, long-sleeved clothing, and suitable splash goggles
--	--

	or face shields when handling batteries; 7. If the slot or cover of the battery is cracked, deformed, or leaks, replace it immediately.
--	---

4.3 System Diagnosis

This product will warn the abnormal state of the vehicle through the flashing of the controller LED light and the horn looping alarm sound. Once the warning signal is generated, it indicates that abnormal condition occurs in some part of the electric wheelchair, which may occur in the motor, brake, battery, wire connection and the product itself. The product has built-in information consultation function, and the displayed diagnosis number can reflect the nature of the condition. These abnormal conditions can still be detected without using other service tools. Voice message indicates: The first two high-frequency sounds are the guiding sounds, and then are the number of alarm sounds. The circular alarm may stop using the electric wheelchair according to the severity of the situation. Please do not continue to use it, and troubleshoot first to ensure that safety.

4.3.1 Sound signal code and fault maintenance light indication: Through the built-in information diagnosis function of the product, the displayed diagnosis number can reflect the nature of the situation;

Code	Problem Description	Diagnosis	Recommended Action
1	All the LED lights are off	Power off or sleep mode	Turn on or check the fuse
2	Fewer LED lights mean less battery power	Low power and low voltage	Low battery voltage, use after charging, if battery damaged replenishing battery damaged, use after replacing. (Please pay attention to environmental protection.
3	The left drive wheel is not moving	Left side motor failure	Check if the left motor or the linking device and motor wires connection are loose.
4	The right drive wheel is not moving	Right side motor failure	Check if the right motor or the linking device and motor wires connection are loose.

5	The controller wire is abnormally hot	Over current	If the motor drive mechanism is stuck, use the ammeter to detect the current and check the controller
6	Direction control is no use	Joystick failure	If the rocker is reset or the rocker cable is disconnected, check whether the link plug is loose.
7	Controller not responding	Controller failure	Please contact with the after-sales service
8	Direction control error	Signal error	Please contact with the after-sales service

4.4 Control System

Please do not use a wheelchair until you are fully familiar with how to use this controller. All controller parameters are professionally set by the factory to produce the best performance under various conditions. If the controller parameters need to be changed, must have the company's professionals to adjust and confirm the user's safety and comply with the safety regulations of medical devices. The circuit design on the wheelchair is according to your safe use. Please do not change or adjust it in any way to cause damage.

4.4.1 Control Panel: 【See Figure 24】




Figure 24

4.4.2 Controller power switch and sleep function

- © Power on the controller: turn on the power switch button  , the power indicator lights up from left to right, and the buzzer beeps, the controller is normally Power on;
- © Power off the controller: Press the power switch button  , all LED lights on the panel are off, and the controller has no output. In certain emergency situations, you can press the power switch button to turn off the controller;
- © Sleep mode: After 20 minutes without any operation of the joystick, the controller automatically turns off the power and enters the sleep mode. Press the key again to wake up the control system. In the sleep state of the product, its own power consumption ($\leq 2.2\text{mA}$) will not damage the battery.

4.4.3 Basic operation of the controller:

4.4.3.1 Make sure the joystick of the controller is perpendicular to the ground, press the power on button on the controller, and the power indicator will display. If the joystick moves first, and then press the power on button, the safety device will make the wheelchair unable to move, and the power display light will be in a flashing state. As long as you release the joystick, it will move;

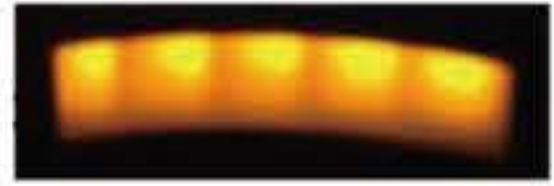4.4.3.2 Speed adjustment: The speedometer displays the comparison value between the current speed and the maximum speed, and the speed can be increased (decreased) by pressing the speed up button  (speed down button ); Each press will increase(decrease) 20%. For fine-tuning, press and hold the button until the speedometer jump to the desired speed, and then release. The speed is divided into five gears  (one light is 20% of the highest speed, five lights for the highest speed);

4.4.3.3 According to the speed limit set by the user, the joystick can control the travel direction [as shown in Figure 25]; the movement direction of the joystick instructs the wheelchair to move, and the moving amplitude of the joystick determines the speed of the wheelchair in this direction:

- Joystick up: drive forward;
- Joystick down: drive backwards;
- Joystick left: turn left;
- Joystick right: turn right.

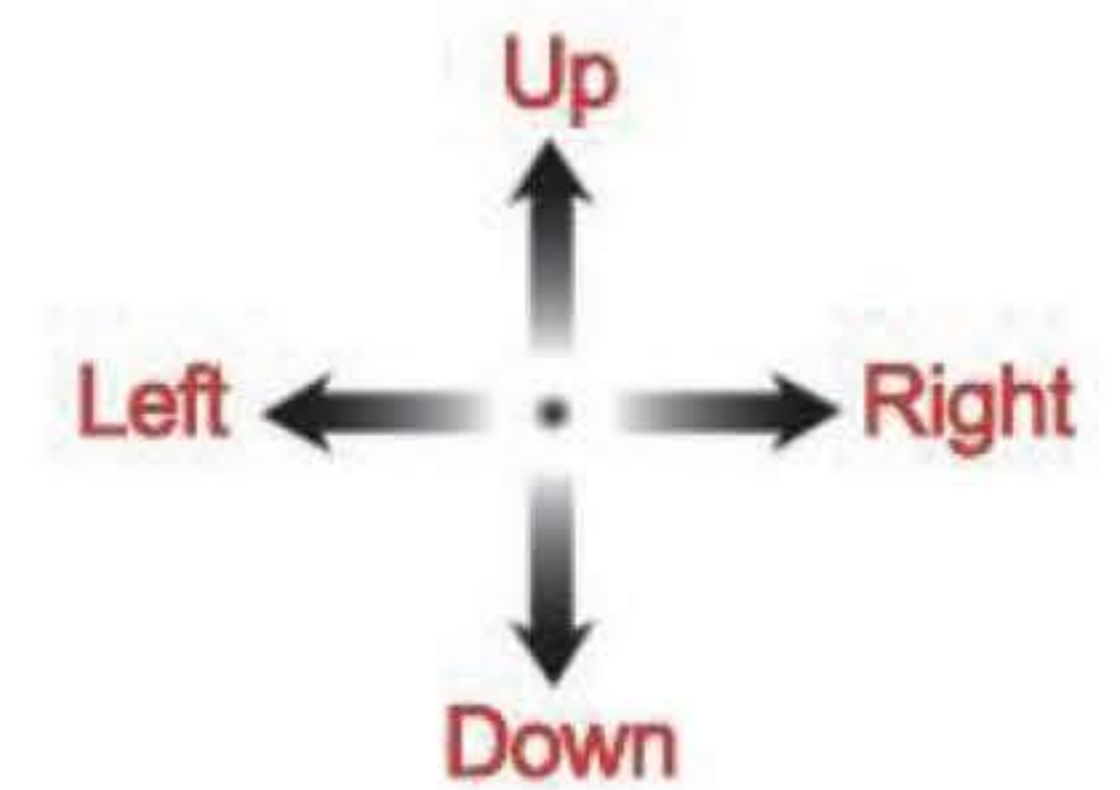


Figure 25

- 4.4.3.4 When start or stop the wheelchair, please adjust the speed adjustment button to a slow speed;
- 4.4.3.5 When driving more proficient, increase the speed to a safe controllable speed;
- 4.4.3.6 For indoor use, please use the lowest speed to avoid collision with objects;
- 4.4.3.7 In case of emergency stop, just release the joystick, or turn it to the opposite direction (emergency stop).

4.4.4 Controller connection diagram

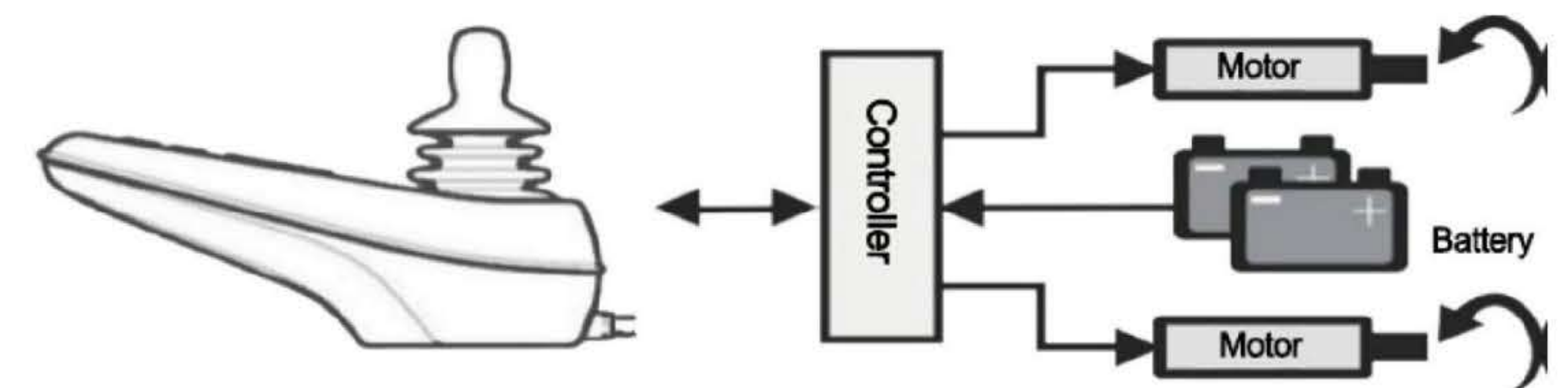


Figure 26

5. Maintenance

5.1. Cleaning and routine inspection:

5.1.1 Clean the wheelchair;

- ◎ Use a wet cloth and mild detergent for cleaning;
- ◎ The electronic part should not be contact with water;
- ◎ Do not use high pressure cleaning equipment;
- ◎ Manual brake: Before each use, make sure that the static function of the wheel after manual brake is normal;
- ◎ Tire pressure: Before each use, be sure to check the tire pressure is normal;
- ◎ Cushions and backrests: Wash with warm water and diluted soapy water if necessary;
- ◎ Controller: In any case, the protective measures of the controller must be taken;
- ◎ Fastening parts: Before each use, check whether all screws are firm, and tighten it immediately if it is loose.

5.1.2 Routine inspection:

Please check the following items before driving, if any abnormality, please contact our after-sales service or the dealer you purchased to get the best solution.

Item	Checking list
Controller	◎ Whether it can be switched on and off normally ◎ Whether the joints or parts are loose ◎ Whether it is easy to operate in any direction ◎ Whether the speed can be adjusted ◎ Whether the power display light is on and there is enough power for use ◎ Whether the horn can be used normally ◎ Is there any abnormal light ◎ Is there any abnormal sound

Motor	◎ Is there any abnormal sound ◎ Whether the electromagnetic brake can work normally (The wheelchair starts and stops with a click).
Clutch	◎ Is there any abnormal sound ◎ Whether the clutch device can work normally
Frame	◎ Is there any abnormal sound ◎ Whether the parts are loose
Tire	◎ Whether the parts are loose ◎ Whether there is damage or insufficient tire pressure ◎ Is there any abnormal sound



Note

If you find any abnormality, please contact our customer service staff or the dealer you purchased for repair service

5.2 Regular maintenance records:

5.2.2 In order to ensure that your wheelchair is in good condition, please contact our customer service staff regularly to do further wheelchair maintenance and make regular maintenance records;

5.2.3 We recommend that your wheelchair be inspected every 6 months and repaired once a year;

5.2.4 The following is the check list for wheelchairs. Please check your wheelchair at the recommended frequency below. For some checks, you must check yourself when you sit in or get up from the wheelchairs to ensure safety.

5.2.5 In order to make you pay more attention, we have listed these self-examination parts as:

A: Check weekly **B:** Check monthly **C:** check every 6 months.
D: Check annually

The following table lists the inspection items. For the health and safety of you and the people around you, please be sure to implement them.

A Check weekly	Check whether the following items are abnormal: ◎ Whether the battery functions normally; ◎ Whether the front and rear wheel parts are normal or whether there is abnormal noise; ◎ Whether the pelvis belt is stable; ◎ Whether the stability of the wheelchair structure is normal; ◎ Whether the activity part is lubricated(Add oil or grease).
B Check monthly	Check the following items for looseness or wear: ◎ Screws and joystick devices of the armrest set; ◎ Electromagnetic brake and manual brake function; ◎ Clutch function; ◎ Joystick; ◎ Front wheel and rear wheel parts; ◎ Depth of front/rear tire tread; ◎ Confirm the wiring of the electronic control system; ◎ The connector between the charger and the controller.
C Check every 6 months.	Please send the wheelchair to the service center for professional inspection once.
D Check annually	It is recommended to send the wheelchair to the service center for professional examination once

6. After-sales service

6.1 Product Term:

The electric wheel chair produced by our company belongs to active, non-sterile, reusable products, and does not contain chemical substances. Through searching, there are no relevant regulations or standards to specify the contents related to the validity period of electric wheel chair and similar products.

The safety, validity period and service life of wheelchair in clinical use depend not only on the structural strength of the product itself, but also on the user's use, the product's use environment, user's use habits, whether it is maintained according to the cycle and other factors. The suggested validity period of our products is 5 years from the date of

production.

6.2 Product Warranty Period:

- 1) Provide 2-year warranty service for the main frame;
- 2) Within 1 year from the date of purchase, if any material or manufacturing defects are found in the parts listed below after inspection by the sales representative, the company will provide free repair or replacement services for the original customers.
 - Electrical control or joystick system;
 - Motor/drive system;
 - Bearings and bushings.
- 3) Parts and conditions not covered by warranty:
 - Consumables are not covered by the warranty, such as: batteries, tires, seat back pads, guard plates, armrest pads and other plastic parts, rubber parts, and other auxiliary equipment;
 - Damage caused by abuse, misuse, accident and negligence;
 - Damage caused by incorrect operation, maintenance and storage;
 - Commercial use or other abnormal use;
 - Vehicle damage caused by various natural disasters and man-made accidents;
 - Sensory judgments that can be identified as having no functional impact;
 - Surface deterioration phenomenon of wear and tear after use;
 - Unauthorized modification of electric wheelchairs;
 - If you do not contact the company or after-sales service outlets, the costs incurred by self-maintenance shall be borne by oneself.
- 4) Comply with maintenance items:
 - According to the operation method shown in this manual, operate normally;
 - Perform the daily inspection and maintenance work shown in this manual.
- 5) Warranty and effective
 - The after-sales service card of the product has been included in this manual, and the after-sales service card will take effect from the date of purchase.
 - The main frame of this product is guaranteed for 2 years (from the date of production), other parts refer to the above regulations.

7. Electromagnetic compatibility

The following cables must be used to comply with electromagnetic emissions and interference resistance requirements

7.1 In addition to the accessories and cables provided by the manufacturer of the electric wheelchair, the use of accessories and cables outside the regulations may lead to an increase in the emission of the electric wheelchair or a decrease in the noise immunity.

7.2 The electric wheelchair should not be used close to or stacked with other equipment. If it must be used close to or stacked, it should be observed and verified that it can operate normally under the configuration it uses.

7.3 The following battery chargers and cables must be used to meet the requirements of electromagnetic emission and anti-interference.

(See Table 1 and Table 2) : Table 1 Charger information

Charger:	24V 2A
----------	--------

Table 2 Cable information

Cable Type	Length
Power cord	0.9m
Charging cord	0.95m

7.4 Except for cables (transducers) sold as spare parts for internal components, the use of unspecified accessories and cables (transducers) may increase the electromagnetic emission or reduce the immunity of the electric wheelchair.

7.5 The electric wheelchair has unique identification of the wheelchair and battery charger model or other information

7.6 Guidelines and manufacturer's declaration.

Table 3 Electromagnetic emission — impact on wheelchair drive system

Launch test	Compliance	Electromagnetic environment — Guidelines
RF emission CISPR 11	Group 1	The equipment only uses RF energy for its internal functions, so its RF emission is very low and there is little possibility of interference to nearby electronic equipment
RF emission CISPR 11	Class B	The equipment is suitable for use in all facilities, including domestic and residential public low-voltage power supply networks directly connected to domestic use.
<u>Harmonic currents emissions</u> <u>IEC 61000-3-2</u>	<u>Class A</u>	<u>The equipment is suitable for use in all facilities, including domestic and residential public low-voltage power supply networks directly connected to domestic use.</u>
Voltage fluctuation /flicker emissions	<u>Compliant</u>	<u>The equipment is suitable for use in all facilities, including domestic and residential public low-voltage power supply networks directly connected to domestic use.</u>

IEC61000-3-3

Table 4 Electromagnetic immunity to wheelchair drive system

Immunity test	Test level	Compliance level	Electromagnetic environment — Guidelines
Electrostatic discharge (ESD) IEC 61000-4-2	ISO 7176-21:2009 $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 6\text{kV}$ contact discharge $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 8\text{kV}$ air discharge IEC 60601-1-2:2014+A1:2020 $\pm 8\text{kV}$ contact discharge $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 8\text{kV}$ $\pm 15\text{kV}$ air discharge	ISO 7176-21:2009 $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 6\text{kV}$ contact discharge $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 8\text{kV}$ air discharge IEC 60601-1-2:2014+A1:2020 $\pm 8\text{kV}$ contact discharge $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 8\text{kV}$ $\pm 15\text{kV}$ air discharge	The floor shall be wood, concrete or ceramic tile, and if the floor is covered with synthetic material, the relative humidity shall be at least 30%
Power frequency magnetic field (50/60Hz) IEC 61000-4-8	30A/m @50Hz or 60Hz	30A/m @50Hz or 60Hz	The power frequency magnetic field shall have the horizontal characteristics as in typical places of commercial or hospital environment

Guidelines and manufacturer's declaration - Electromagnetic Immunity

The purchaser or user shall ensure that the equipment is intended for use in the electromagnetic environment specified below:

Immunity test	Test level	Coincidence level	Electromagnetic environment — Guidelines
Radio frequency radiation IEC 61000-4-3	ISO 7176-21:2009 20V/m @26MHz-2,5GHz AM 80% at 1 kHz IEC 60601-1-2:2014+A1:2020 3V/m @80MHz-2,7GHz AM 80% at 1 kHz	ISO 7176-21:2009 20V/m @26MHz-2,5GHz AM 80% at 1 kHz IEC 60601-1-2:2014+A1:2020 3V/m @80MHz-2,7GHz AM 80% at 1 kHz	Portable and mobile RF communication equipment should not be used closer to any part of the equipment than the recommended isolation distance, including cables. The distance is calculated by the formula corresponding to the transmitter frequency. Recommended isolation distanced $d = 1.2\sqrt{(P)}$ 150kHz-80MHz $d = 1.2\sqrt{(P)}$ 80MHz-800MHz $d = 2.3\sqrt{(P)}$ 800MHz-2.5GHz Note : P - according to the maximum rated output power of

			the transmitter provided by the transmitter manufacturer, in watts (W); D - is the recommended isolation distance, in meters (m). The field strength of the fixed RF transmitter is determined by surveying the electromagnetic field A, and in each frequency range, it shall be lower than the coincidence level at B. Interference may occur near equipment marked with the following symbols. 
A Stationary transmitter, such as base stations of wireless (cellular / cordless) telephones and ground mobile radios, amateur radios, AM and FM radio broadcasting and television broadcasting, have field strengths that cannot be predicted accurately in theory. In order to evaluate the electromagnetic environment of fixed RF transmitter, the survey of electromagnetic field should be considered. If the measured field strength of the place where the equipment is located is higher than the above applicable RF compliance level, the equipment shall be observed to verify its normal operation. If abnormal performance is observed, supplementary measures may be necessary, such as reorienting or re-positioning the equipment. B in the whole frequency range of 150kHz-80mhz, the field strength shall be less than 3V / m.			

Recommended isolation distance between portable and mobile RF communication equipment. The equipment is expected to be used in an electromagnetic environment with controlled RF radiation disturbance. According to the maximum rated output power of communication equipment, the buyer or user can prevent electromagnetic interference by maintaining the minimum distance between portable and mobile RF communication equipment (transmitter) and equipment as recommended below:			
Maximum rated output power of transmitter W	The isolation distance corresponding to different frequencies of the transmitter/m		
	150kHz - 80MHz $d = 1.2\sqrt{(P)}$	80MHz - 800MHz $d = 1.2\sqrt{(P)}$	800MHz - 2.5GHz $d = 2.3\sqrt{(P)}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For the maximum rated output power of the transmitter not listed in the above table, the recommended isolation distance D, in meters (m), can be determined by the formula in the corresponding transmitter frequency column, where p is the maximum rated output power of the transmitter provided by the transmitter manufacturer. In watts (W). Note 1: At 80MHz and 800MHz frequency points, the formula of higher frequency band is adopted. Note 2: These guidelines may not be suitable for all situations. Electromagnetic transmission is affected by the absorption and reflection of buildings, objects and human bodies.			

Table 5 Electromagnetic emission - for vehicle mounted and on-vehicle mounted battery chargers		
Guidelines and manufacturer's statements - electromagnetic emissions		
The purchaser or user shall ensure that the equipment is intended for use in the electromagnetic environment specified below:		
Launch test	Compliance	Electromagnetic environment—Guidelines
RF emission CISPR 11	Group 1	The equipment only uses RF energy for its internal functions, so its RF emission is very low and there is little possibility of interference to nearby electronic equipment
RF emission CISPR 11	Class B	The equipment is suitable for use in all facilities, including domestic and residential public low-voltage power supply networks directly connected to domestic use.
Harmonic radiation IEC 61000-3-2	Class A	
Voltage fluctuation /flicker emission IEC 61000-3-3	Compliance	

Table 6 Electromagnetic immunity - for vehicle mounted and non-vehicle mounted battery chargers			
Guidelines and manufacturer's statements - electromagnetic immunity			
The purchaser or user shall ensure that the equipment is intended for use in the electromagnetic environment specified below:			
Immunity test	Test level	Compliance level	Electromagnetic environment—Guidelines
Electrostatic discharge (ESD) IEC 61000-4-2	ISO 7176-21:2009 $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 6\text{kV}$ contact discharge $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 8\text{kV}$ air discharge IEC 60601-1-2:2014+A1:2020 $\pm 8\text{kV}$ contact discharge $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 8\text{kV}$ $\pm 15\text{kV}$ air discharge	ISO 7176-21:2009 $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 6\text{kV}$ contact discharge $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 8\text{kV}$ air discharge IEC 60601-1-2:2014+A1:2020 $\pm 8\text{kV}$ contact discharge $\pm 2\text{kV}$ $\pm 4\text{kV}$ $\pm 8\text{kV}$ $\pm 15\text{kV}$ air discharge	Floors should be wood, concrete or tile, if floors are covered with synthetic material, the relative humidity should be at least 30%
Group of electrical fast transient pulses IEC 61000-4-4	ISO 7176-21:2009 $\pm 1\text{kV}$ Power cord $\pm 0.5\text{kV}$ for input/output lines IEC 60601-1-2:2014+A1:2020 $\pm 2\text{kV}$ Power cord $\pm 1\text{kV}$ for input/output lines 100 kHz repetition frequency	ISO 7176-21:2009 $\pm 1\text{kV}$ Power cord IEC 60601-1-2:2014+A1:2020 $\pm 2\text{kV}$ Power cord 100 kHz repetition frequency	Network power supply should be of typical quality used in a commercial or hospital environment.

Surge IEC 61000-4-5	<u>ISO 7176-21:2009</u> ± 1kV Cord to cord ± 2kV cord to the ground <u>IEC 60601-1-2:2014+A1:2020</u> ± 0.5kV ± 1kV Cord to cord ± 0.5kV ± 1kV ± 2kV cord to the ground	<u>ISO 7176-21:2009</u> ± 1kV Cord to cord ± 2kV cord to the ground <u>IEC 60601-1-2:2014+A1:2020</u> ± 0.5kV ± 1kV Cord to cord ± 0.5kV ± 1kV ± 2kV cord to the ground	Network power supply should be of typical quality used in a commercial or hospital environment.
Voltage dips, short interruptions, and voltage variations on power input lines IEC 61000-4-11	<u>ISO 7176-21:2009</u> 100% UT, Lasting 0.5 cycle, 100% UT, lasting one cycle 30% UT, lasting 25 cycle <u>IEC 60601-1-2:2014+A1:2020</u> 0 % UT; 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0° 0 % UT; 250/300 cycle	<u>ISO 7176-21:2009</u> 100% UT, Lasting 0.5 cycle, 100% UT, lasting one cycle 30% UT, lasting 25 cycle <u>IEC 60601-1-2:2014+A1:2020</u> 0 % UT; 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0° 0 % UT; 250/300 cycle	The network power supply shall be of the quality typical for use in a commercial or hospital environment. If the user needs to run continuously during the power interruption, it is recommended that the device be powered by an uninterrupted power supply or battery

Note: UT refers to the AC network voltage before applying the test voltage.

Guidelines and manufacturer's statements - electromagnetic immunity			
The purchaser or user shall ensure that the equipment is intended for use in the electromagnetic environment specified below:			
Immunity test	Test level	Compliance level	Electromagnetic environment—Guidelines
RF transmission IEC 61000-4-6 RF radiation IEC 61000-4-3	3V (effective value) 150kHz - 80MHz 3V/m 80MHz - 2.5GHz	3V (effective value) 150kHz - 80MHz 3V/m 80MHz - 2.5GHz	Portable and mobile RF communications equipment should not be used closer than the recommended isolation distance to any part of the equipment, including cables. The distance is calculated by the formula corresponding to the frequency of the transmitter Recommended isolation distance

			$d = 1.2\sqrt{P}$ 150kHz-80MHz $d = 1.2\sqrt{P}$ 80MHz-800MHz $d = 2.3\sqrt{P}$ 800MHz-2.5GHz Note: P - according to the maximum rated output power of the transmitter provided by the transmitter manufacturer, in watts (W); D - is the recommended isolation distance, in meters (m). The field strength of the fixed RF transmitter is determined by surveying the electromagnetic field A, and in each frequency range, it shall be lower than the coincidence level at B. Interference may occur near equipment marked with the following symbols. 
A Stationary transmitter, such as base stations of wireless (cellular / cordless) telephones and ground mobile radios, amateur radios, AM and FM radio broadcasting and television broadcasting, have field strengths that cannot be predicted accurately in theory. In order to evaluate the electromagnetic environment of fixed RF transmitter, the survey of electromagnetic field should be considered. If the measured field strength of the place where the equipment is located is higher than the above applicable RF compliance level, the equipment shall be observed to verify its normal operation. If abnormal performance is observed, supplementary measures may be necessary, such as reorienting or re-positioning the equipment. B in the whole frequency range of 150kHz-80mhz, the field strength shall be less than 3V / m.			

Recommended isolation distance between portable and mobile RF communication equipment. The equipment is expected to be used in an electromagnetic environment with controlled RF radiation disturbance. According to the maximum rated output power of communication equipment, the buyer or user can prevent electromagnetic interference by maintaining the minimum distance between portable and mobile RF communication equipment (transmitter) and equipment as recommended below:			
Maximum rated output power of transmitter W	The isolation distance corresponding to different frequencies of the transmitter/m		
	150kHz - 80MHz $d = 1.2\sqrt{P}$	80MHz - 800MHz $d = 1.2\sqrt{P}$	800MHz - 2.5GHz $d = 2.3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For the maximum rated output power of the transmitter not listed in the above table, the recommended isolation distance D, in meters (m), can be determined by the formula in the corresponding transmitter frequency column, where p is the maximum rated output power of the transmitter provided by the transmitter manufacturer. In watts (W).

Note 1: At 80MHz and 800MHz frequency points, the formula of higher frequency band is adopted.

Note 2: These guidelines may not be suitable for all situations. Electromagnetic transmission is affected by the absorption and reflection of buildings, objects and human bodies.

Guidelines and manufacturer's statements - electromagnetic immunity

The purchaser or user shall ensure that the equipment is intended for use in the electromagnetic environment specified below:

Immunity test	Test level	Compliance level	Electromagnetic environment— —Guidelines
Proximity fields from RF wireless communications Equipment IEC 61000-4-3	3V/m @380-390MHz Pulse Modulation 18Hz 28V/m @430-470MHz FM ± 5 kHz deviation 1 kHz sine 9V/m @704-787MHz Pulse Modulation 217Hz 28V/m @800-960MHz Pulse Modulation 18Hz 28V/m @1700-1990MHz Pulse Modulation 217Hz 28V/m @2400-2570MHz Pulse Modulation 217Hz 9V/m @5100-5800MHz Pulse Modulation 217Hz	3V/m @380-390MHz Pulse Modulation 18Hz 28V/m @430-470MHz FM ± 5 kHz deviation 1 kHz sine 9V/m @704-787MHz Pulse Modulation 217Hz 28V/m @800-960MHz Pulse Modulation 18Hz 28V/m @1700-1990MHz Pulse Modulation 217Hz 28V/m @2400-2570MHz Pulse Modulation 217Hz 9V/m @5100-5800MHz Pulse Modulation 217Hz	The distance between the microwave therapeutic system and RF wireless communications equipment shall be greater than or equal to 30 cm (12 inches).

WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the microwave therapeutic system, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

Attention !

Our wheelchairs and battery chargers meet the requirements of GB/TISO 7176-21:2009 "Electromagnetic compatibility Requirements and Test Methods for Electric wheelchairs, Electric mobility Scooters and Battery Chargers" and IEC 60601-1-2:2014+A1:2020 "Medical Electrical Equipment -- Part 1-2: General Safety Requirements and Standards: Electromagnetic Compatibility Requirements and Tests".

8. Packing list

NO.	Item name	Quantity	Note
1	Wheelchair (with foot pedal)	1unit	Mainframe
2	Battery pack	1pc	Removable
3	Charger	1pc	/
4	Instruction manual (including after-sales service card)	1pc	/
5	Upper seat cushion, upper back cushion, leg straps	1set	Detachable

9. Power supply conditions and environmental requirements

Working environment requirements

Temperature 5 °C ~40 °C , relative humidity 15%~90% and atmospheric pressure range 700hPa~1060hPa。

Working power supply conditions

Charging power supply: 220 V 50 Hz, maximum output current 2 A.



Internal power supply device operating voltage: DC 24 V.

Transportation and storage conditions:

Storage and transportation temperature should not exceed -20℃~55℃, humidity should not exceed 93%, atmospheric pressure range 700hPa~1060hPa, no corrosive gas and clean room with good ventilation.

Loading and transporting according to the transportation marks and graphics on the outsourcing box.

10.Labels,Packing Logo Design

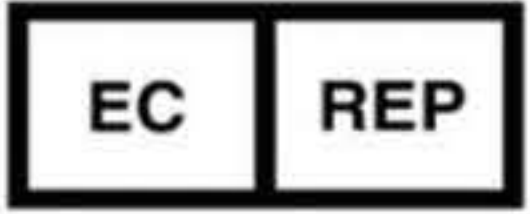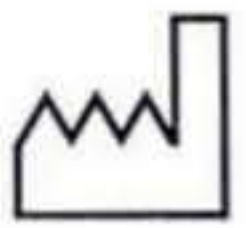Symbol	Introductions	Symbol	Introductions
	CE Symbol		AUTHORISED REPRESENTATIVE IN THE EUROPEAN COMMUNITY
	Manufacturer Name Address		Manufacture Date
	Batch Code		Store in clean & dry place protected from rain, snow, ice, salt and water.
	Warnings and Precautions		non-sterile
	medical device		Use-by date
	Keep dry		Keep away from Sunlight
	Consult instructions for use		Unique Device Identification
	Country of Manufacture		Model Number
	Importer		Distributor